

Maximizer[™]24 Green

Single Use Gamma Irradiated Hollow Fiber Cartridges. Easy and Reliable



The Maximizer[™] Green production scale hollow fiber cartridges are designed for biopharmaceutical concentration/diafiltration and clarification unit operations. With a surface area of 5.42 m², the Maximizer[™]24 Green HF cartridge can accommodate batch volumes between 100L and 500L in a single cartridge system, and several thousand liters with a multi-cartridge set up. It has an average process flow rate of 153-306 L/hr for ultrafiltration applications and 102-204 L/hr for microfiltration applications.


The Maximizer[™] Green cartridges are gamma irradiated and ready to use, and utilize WaterSep's antifouling, low binding, modified polyethersulfone membrane (m-PES) which provides a process flux and product recovery that typically exceeds other crossflow devices.

The Maximizer[™]24 Green hollow fiber cartridges are offered with the same path length as other high performance WaterSep HF cartridges, which makes scale up and scale down easy and predictable.

Maximizer[™] Green cartridges are also offered as part of fully integrated open architecture single use crossflow assemblies available with pressure sensors, flow meters, pumpheads, aseptic connectors or bags.

The Maximizer[™]24 Green cartridges offer a true single use crossflow system that provides:

- Significant cost savings over traditional cassette formats.
- Gamma irradiated - ready to use - no pre-rinse required.
- No humectants - 80x less extractables, see [Extractables Flush Analysis](#).
- No need for time-consuming cleaning validation.
- Reduced risk of cross contamination.
- Consistent membrane performance batch to batch.
- Self-contained – no assembly – no expensive hardware – no membrane installation.
- Every HF cartridge is integrity tested and has an individual lot number for easy traceability.

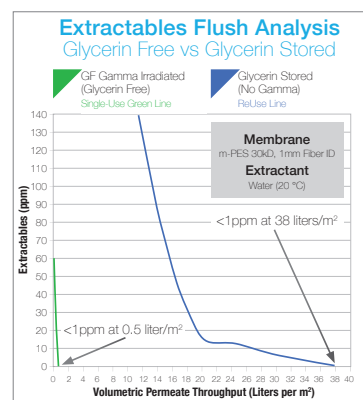

TYPICAL APPLICATIONS

- Clarification of **mammalian/CHO** cell cultures and maximizing **protein recovery**.
- **Cell-harvest**. (Excellent results have been achieved with both *E. Coli* whole cells and *E. Coli* lysates, as well as other microbial process streams.)
- Concentration and purification of **vaccines**.
- Concentration and diafiltration of **gene therapy** products.
- Concentration/diafiltration of **monoclonal antibodies, recombinant proteins, biological macromolecules and peptides**.

Maximizer™ 24 Green

ORDERING INFORMATION

Place orders online at www.watersep.net, or email purchase orders to orders@watersep.net, sales order confirmation to follow.

Maximizer™ 24 Green HF Cartridge

1.0 mm Fiber ID
SU XXX 10MAX24SH, 1/pkg

Maximizer™ 24 Green HF Cartridge

2.0 mm Fiber ID*
SU XXX 20MAX24 SH, 1/pkg
*Available in 10K, 50K, 300K-750K and 0.1 μm - 0.45 μm

Replace XXX with:

003 for 3K	050 for 50K	750 for 750K
005 for 5K	100 for 100K	910 for 0.1 μm
010 for 10K	300 for 300K	920 for 0.2 μm
030 for 30K	500 for 500K	945 for 0.45 μm

SPECIFICATIONS

Cartridge Dimensions
4.62" (11.73 cm) x 28.80" (73.15 cm)

Membrane Surface Area
58.4 ft² (5.42 m²)

Molecular Weight Cut-off
3K, 5K, 10K, 30K, 50K, 100K, 300K, 500K, 750K

Membrane Pore Size
0.1 μm , 0.2 μm and 0.45 μm

Fiber ID
1.0 mm and 2.0 mm

MATERIALS OF CONSTRUCTION

The cartridge as assembled is USP Class VI compliant.

Membrane
Modified Polyethersulfone (m-PES)


Housing
White Polysulfone

Encapsulant
Medical Grade Epoxy

CONNECTIONS

Feed/Retentate
2.0" TC Standard Sanitary Fitting

Permeate
1.5" TC Standard Sanitary Fitting

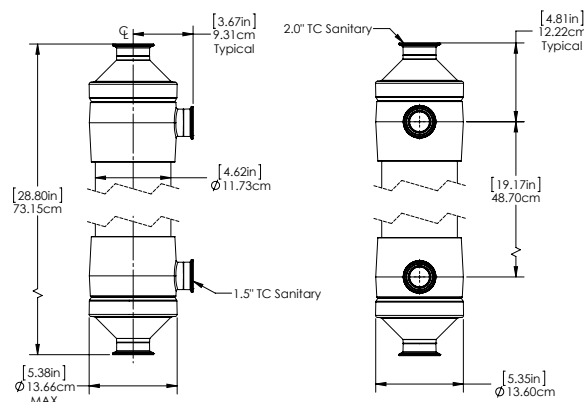

Technical Support: experttalk@watersep.net
or call 508.970.0089 x 204.