

## Maximizer™ 41 Green

### Single Use Gamma Irradiated Hollow Fiber Cartridges. Easy and Reliable



The Maximizer™ Green production scale hollow fiber cartridges are designed for biopharmaceutical concentration/diafiltration and clarification unit operations. With a surface area of 9.74 m<sup>2</sup>, the Maximizer™ 41 Green HF cartridge can accommodate batch volumes between 300L and 1,000L in a single cartridge system, and several thousand liters with a multi cartridge set up. It has an average process flow rate of 240-480 L/hr for ultrafiltration applications and 160-320 L/hr for microfiltration applications.


The Maximizer™ Green cartridges are gamma irradiated and ready to use, and utilize WaterSep's antifouling, low binding, modified polyethersulfone membrane (m-PES) which provides a process flux and product recovery that exceeds other crossflow devices.

The Maximizer™ 41 Green hollow fiber cartridges are offered with the same path length as other high performance WaterSep HF cartridges, which makes scale up and scale down easy and predictable.

Maximizer™ Green cartridges are also offered as part of fully integrated open architecture single use crossflow assemblies available with pressure sensors, flow meters, pumpheads, aseptic connectors or bags.

The Maximizer™ 41 Green cartridges offer a true single use crossflow system that provides:

- Significant cost savings over traditional cassette formats.
- Gamma irradiated - ready to use - no pre-rinse required.
- No humectants - 80x less extractables, see [Extractables Flush Analysis](#).
- No need for time-consuming cleaning validation.
- Reduced risk of cross contamination.
- Consistent membrane performance batch to batch.
- Self-contained – no assembly – no expensive hardware – no membrane installation.
- Every HF cartridge is integrity tested and has an individual lot number for easy traceability.

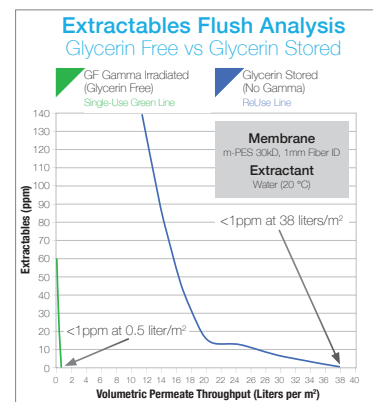

### TYPICAL APPLICATIONS

- Clarification of **mammalian/CHO** cell cultures and maximizing **protein recovery**.
- **Cell-harvest**. (Excellent results have been achieved with both *E. Coli* whole cells and *E. Coli* lysates, as well as other microbial process streams.)
- Concentration and purification of **vaccines**.
- Concentration and diafiltration of **gene therapy** products.
- Concentration/diafiltration of **monoclonal antibodies, recombinant proteins, biological macromolecules and peptides**

## Maximizer™ 41 Green

### ORDERING INFORMATION

Place orders online at [www.watersep.net](http://www.watersep.net), or email purchase orders to [orders@watersep.net](mailto:orders@watersep.net), sales order confirmation to follow.

#### Maximizer™ 41 Green HF Cartridge

1.0 mm Fiber ID  
SU XXX 10MAX41SH, 1/pkg

#### Maximizer™ 41 Green HF Cartridge

2.0 mm Fiber ID\*  
SU XXX 20MAX41 SH, 1/pkg  
\*Available in 10K, 50K, 300K-750K and 0.1  $\mu\text{m}$  - 0.45  $\mu\text{m}$

Replace XXX with:

003 for 3K	050 for 50K	750 for 750K
005 for 5K	100 for 100K	910 for 0.1 $\mu\text{m}$
010 for 10K	300 for 300K	920 for 0.2 $\mu\text{m}$
030 for 30K	500 for 500K	945 for 0.45 $\mu\text{m}$

### SPECIFICATIONS

#### Cartridge Dimensions

4.62" (11.73 cm) x 46.50" (118.11 cm)

#### Membrane Surface Area

104.9 ft<sup>2</sup> (9.74 m<sup>2</sup>)

#### Molecular Weight Cut-off

3K, 5K, 10K, 30K, 50K, 100K, 300K, 500K, 750K

#### Membrane Pore Size

0.1  $\mu\text{m}$ , 0.2  $\mu\text{m}$  and 0.45  $\mu\text{m}$

#### Fiber ID

1.0 mm and 2.0 mm

### MATERIALS OF CONSTRUCTION

The cartridge as assembled is USP Class VI compliant.

#### Membrane

Modified Polyethersulfone (m-PES)

#### Housing

White Polysulfone

#### Encapsulant

Medical Grade Epoxy


### CONNECTIONS

#### Feed/Retentate

2.0" TC Standard Sanitary Fitting

#### Permeate

1.5" TC Standard Sanitary Fitting

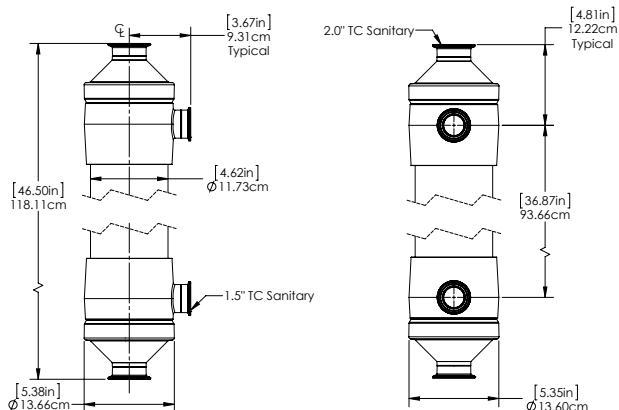

Technical Support: [experttalk@watersep.net](mailto:experttalk@watersep.net)  
or call 508.970.0089 x 204.